
### **i** Opis ogólny / Basic description

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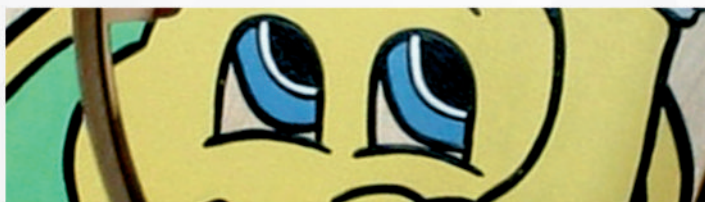
„ELEPHANT” STANDING FRAME is designed for an active rehabilitation of disabled children. The standing frame allows for independent (but under supervision) and safe standing in the appliance, fun and education. Thanks to its simplicity, the „ELEPHANT” STANDING FRAME is an extremely precious device at the rehabilitation of disabled children at home, in hospitals and rehabilitation centers. The correct use of this appliance can help to improve the quality of children’ life.

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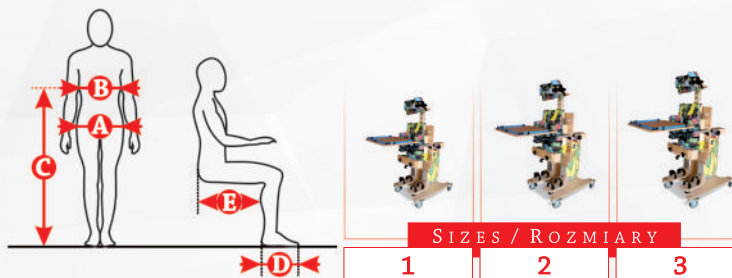	SIZES / ROZMIARY		
	1	2	3
A. Width of hips [cm]	15 - 22	17 - 28	27 - 38
B. Width of the chest [cm]	15 - 22	17 - 28	27 - 38
C. Height from feet to armpit [cm]	48 - 72	58 - 88	68 - 98
D. Foot length [cm][max]	17	20	24
E. Distance from the back to the calf [cm]	11 - 23	19 - 35	24 - 43
1. User weight [kg][max]	25	45	60
2. Dimension of the device [cm][LxWxH]	80x60x65	86x65x75	92x66x80
3. Machine weight [kg]	25	30	35

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## ADJUSTMENT OF HIP CUSHION

In order to adjust the width of cushions supporting the chest, you need to loosen the knob, move the cushion to the desired position and then tighten the knob.

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## CHEST BELT AND HIP BELT

Thanks to the use of strips, belts are easy to fasten.

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## MOVING

The device is equipped with 4 pieces of wheelset. Wheelsets have a ratchet brake.

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## SAFETY BELTS

Adjustment of the safety belt's length is realized by moving the tape in the buckle.

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## CUSHIONS SUPPORTING CHEST AND HIPS


In order to regulate the width of cushions supporting the chest, you need to loosen knobs (located behind the standing frame), adjust cushions to the child's chest and tighten knobs. The position of hip cushions is similarly adjusted.

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## LOWER SUPPORT

This support is used to give the child a full support.

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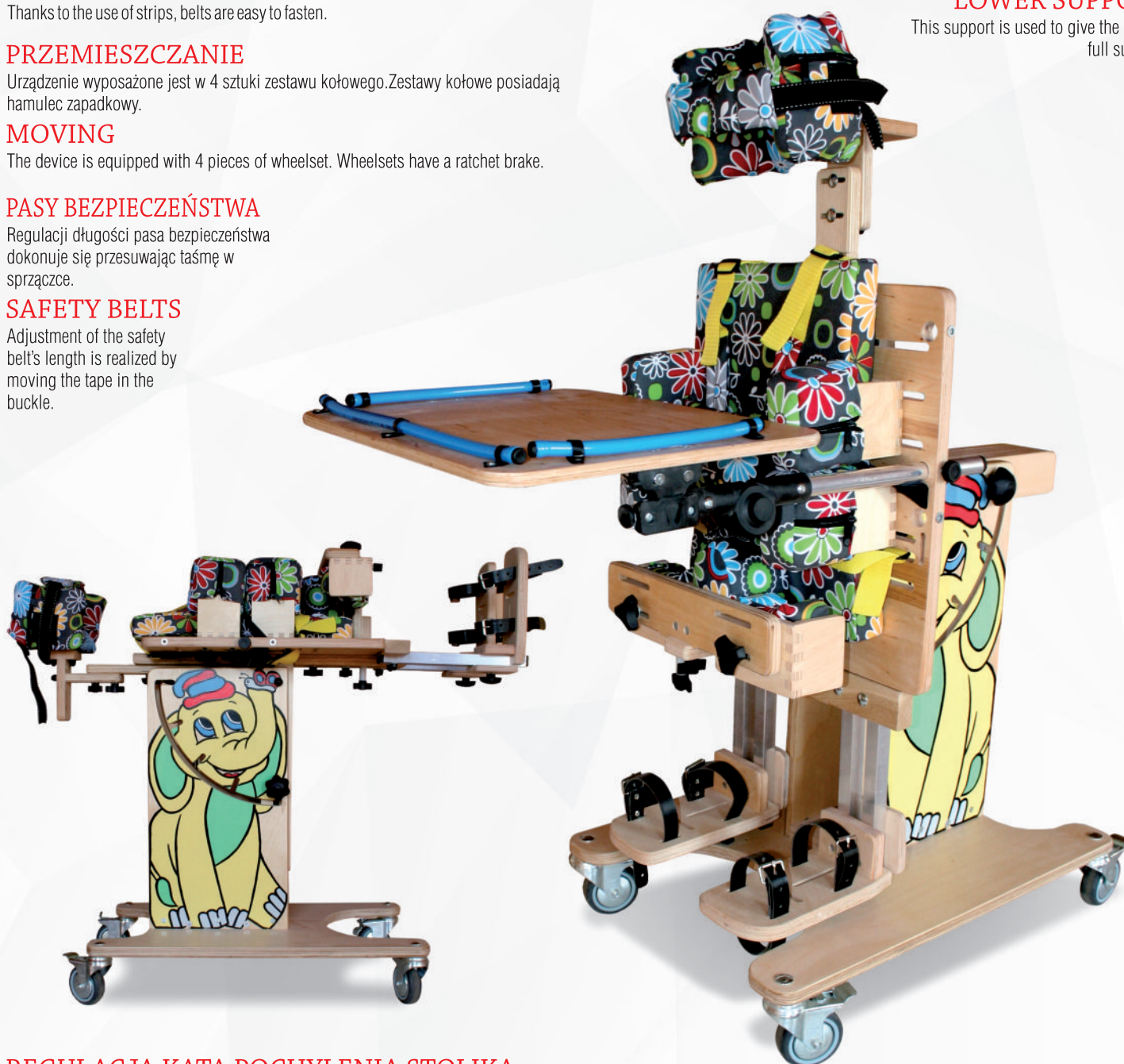## ADJUSTMENT OF TABLE'S INCLINATION ANGLE

By loosening the nuts you can change the table's inclination angle. After setting a desired position, you need to tighten the nuts.

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## VARIATION OF THE STANDING FRAME

In the electric version, the standing frame is adjusted by pressing an appropriate button on the remote control, and in the manual version by turning a proper knob, setting the standing frame and tightening screws.

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## FOOTERS

The position of footers on the footrest plate is adjusted with the use of knobs.